WWW.CEKA-PRECILINE.COM

SWISS MADE

CE 0510



CATALOGUE

CEKA REVAX
CEKA
PRECI-VERTIX
PRECI-VERTIX AT
PRECI-SAGIX
PRECI-BALL
PRECI-CLIX RADICULAR
PRECI-CLIX RADICULAR RC
PRECI-CLIX AXIAL
PRECI-CLIX EXTRACORONAL
PRECI-HORIX
PRECI-BAR
PRECI-CLIP
PRECI-TUBIX/MORTIX
PRECI-POST
3C-BOND
PERMA-RET
EXPANDO
PRECI-PROFILE
CEKA SOL

Prospectus nº 269.E-ed.D
Copyright© ALPHADENT NV 9/2010
BTW BE 0414.789.321 RPR Kortrijk

ARSEUS DENTAL
SERVING THE DENTAL PROFESSIONALS

ALPHADENT NV
Textielstraat 24
8790 Waregem, Belgium
tel +32 56 629 100
fax +32 56 629 101
e-mail: info@ceka-preciline.com

CEKA
ATTACHMENTS
PRECI-LINE


Your partner for precision attachments



CATALOGUE

More at:

WWW.CEKA-PRECILINE.COM

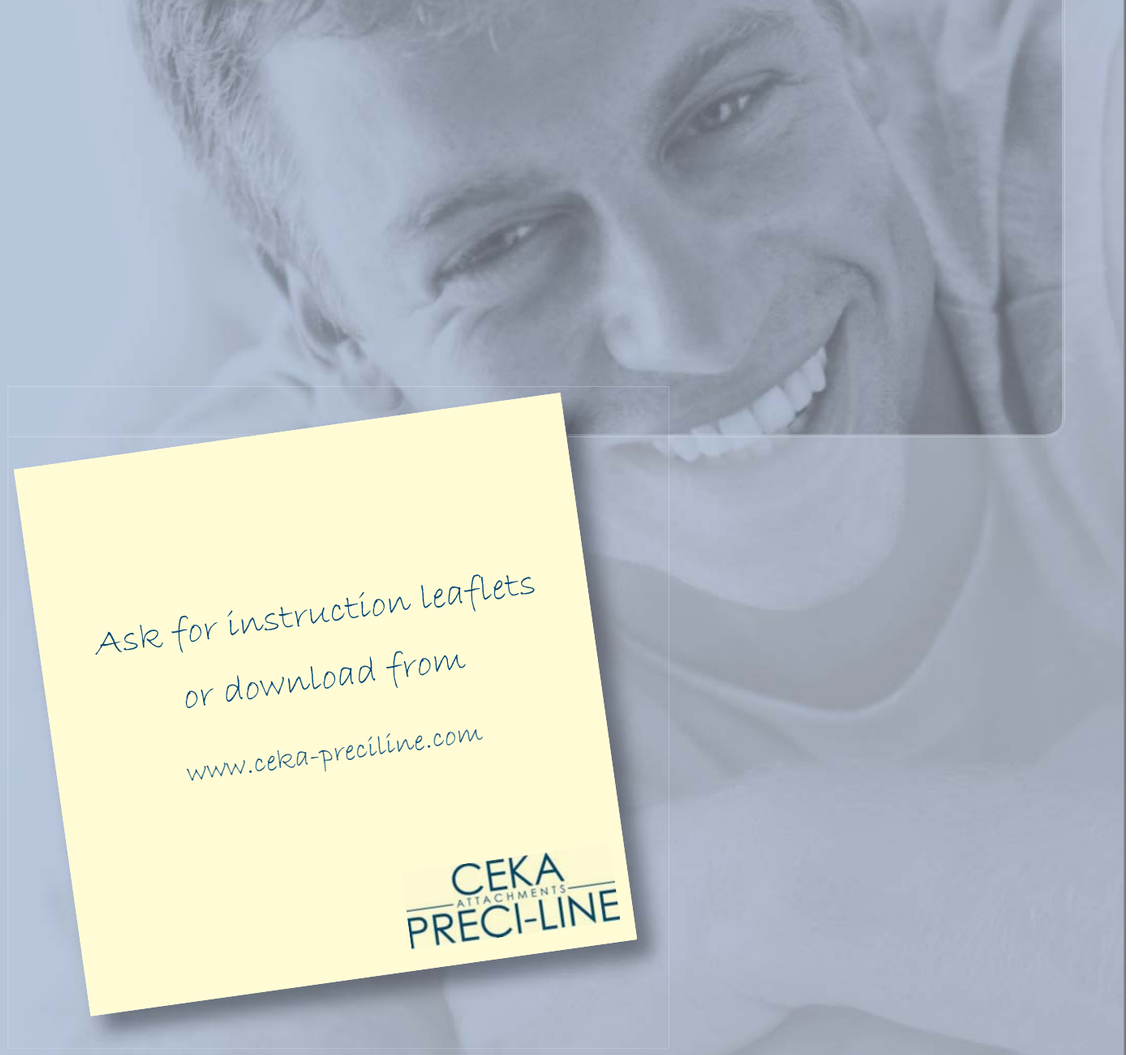

WWW.CEKA-PRECILINE.COM



- Complete catalogue
- All brochures in pdf
- Service and Aftercare
- Safety Data Sheets
- Instructions for use
- FAQ's
- Old spring pins
- ISO/CE certificates
- Promotions

Copyright© ALPHADENT NV 9/2010. No part of this brochure may be reproduced in any manner without permission of ALPHADENT NV. Our products are subject to continuous development. We therefore reserve the right to make technical changes without preliminary notice.

CEKA, CEKA REVAX, PRECI-LINE, PRECI-VERTIX and PRECI-HORIX are registered trademarks of ALPHADENT NV.

CEKA® ATTACHMENTS

M2 size (2 mm)



M2 size male spring pin
 # RE 0031

Height: 2.7 mm
2 pcs + 2 space maintainers

M3 size (3 mm)



M3 size male spring pin
 # 694 C

Height: 3.15 mm
2 pcs + 2 space maintainers


CEKA REVAX – EXTRACORONAL GREEN

M2: H 3.80 mm - ø 3.4 mm
M3: H 4.35 mm - ø 4.4 mm


Female for bonding with CEKA SITE. Use any dental alloy.
Male for bonding with CEKA SITE **(A)**, winged acrylic fixation **(B)**, soldering with CEKA SOL **(C)** or circular acrylic fixation **(D)**.
Each package contains 1 attachment and 3 different plastic profiles (30°, 45°, 60°) to choose from.
See page 6 for individual components, including laser components.



M2 size: **RE 0185 TI**
M3 size: **OL 0285 TI**

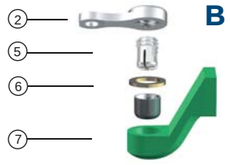

M2 size: **RE 0175 TI**
M3 size: **OL 0275 TI**

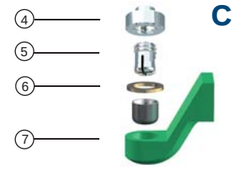

M2 size: **RE 0161 TI**
M3 size: **OL 0261 TI**



M2 size: **RE 0195 TI**
M3 size: **OL 0295 TI**
NEW

CEKA REVAX – EXTRACORONAL BLUE

M2: H 3.80 mm - ø 3.4 mm
M3: H 4.35 mm - ø 4.4 mm

Female for the cast-on technique with non-precious alloys.
Male for bonding with CEKA SITE **(A)**, winged acrylic fixation **(B)**, soldering with CEKA SOL **(C)** or circular acrylic fixation **(D)**.
Each package contains 1 attachment and 3 different plastic profiles (30°, 45°, 60°) to choose from.
See page 6 for individual components, including laser components.



M2 size: **RE 0185 NP**
M3 size: **OL 0285 NP**



M2 size: **RE 0175 NP**
M3 size: **OL 0275 NP**



M2 size: **RE 0161 NP**
M3 size: **OL 0261 NP**



M2 size: **RE 0195 NP**
M3 size: **OL 0295 NP**
NEW

CEKA REVAX – EXTRACORONAL ORANGE

M2: H 3.80 mm - ø 3.4 mm
M3: H 4.35 mm - ø 4.4 mm

Female for the cast-on technique with precious alloys.
Male for bonding with CEKA SITE **(A)**, acrylic fixation **(B)**, soldering with CEKA SOL **(C)** or circular acrylic fixation **(D)**.
Each package contains 1 attachment and 3 different plastic profiles (30°, 45°, 60°) to choose from.
See page 6 for individual components, including laser components.



M2 size: **RE 0185 IR**
M3 size: **OL 0285 IR**



M2 size: **RE 0175 IR**
M3 size: **OL 0275 IR**



M2 size: **RE 0161 IR**
M3 size: **OL 0261 IR**




M2 size: **RE 0195 IR**
M3 size: **OL 0295 IR**
NEW

MORE AT: WWW.CEKA-PRECILINE.COM

CEKA – AXIAL

M2: H 3.6 mm - ø 3.4 mm
M3: H 4.0 mm - ø 4.2 mm

Female for acrylic fixation in the removable prosthesis.
Male on post-coping: for soldering with CEKA SOL FILIGRAN **(A)** or for the cast-on technique with precious alloys **(B)**.
Each package contains 1 attachment. See page 6 for individual components, including laser components.

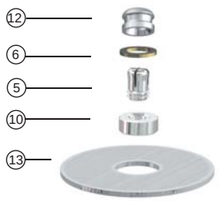

A

M2 size: **RA 61 TI**
M3 size: **691 TI**




B

M2 size: **RA 63 TI**
M3 size: **693 TI**

CEKA – AXIAL & BAR CONSTRUCTIONS


M2: H 3.80 mm - ø 3.4 mm
M3: H 4.35 mm - ø 4.4 mm

Female for bonding with CEKA SITE in post-copings or in bar constructions (all alloys).
Male for bonding with CEKA SITE **(A)** or for acrylic fixation **(B)**.
Each package contains 1 attachment. See page 6 for individual components, including laser components.

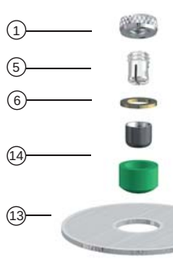

A

M2 size: **RE 0785 TI**
M3 size: **OL 0885 TI**

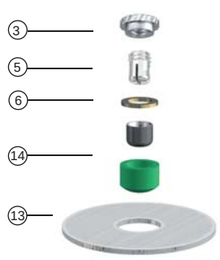


B

M2 size: **RE 0795 TI**
M3 size: **OL 0895 TI**

CEKA ID

M2: H 4.65 mm - ø 3.2 mm
M3: H 4.80 mm - ø 4.0 mm

Intracoronal CEKA combination attachment for the cast-on technique with precious alloys.
A 0.3 mm space maintainer prevents uncontrolled forces from acting on the extracoronal female.
Each package contains 1 attachment. See page 6 for individual components.
Use the P 9 paralleling mandrel.

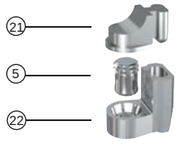

M2 size: **RE 0703 IR**
M3 size: **703 IR**

CEKA 334

H 4.4 mm - ø 4.3 mm

Female for the cast-on technique with non-precious alloys.
Male for the bolt-on technique in the cast metal partial frame.
See page 6 for individual components.



M3 size only: **OL 334 NX**

2 attachments/package



M3 size only: **OL 0234 NPR**

1 attachment/package
Includes 60° profile only.

CEKA CLASSIC

H 3.55 mm - ø 4.3 mm

Female for the cast-on technique with non-precious alloys (blue), precious alloys (orange), or any dental alloy (yellow).
Male for bonding with CEKA SITE **(A)**, acrylic fixation **(B)** or for soldering with CEKA SOL FILIGRAN **(C)**.
Each package contains 2 attachments in one bar profile.
See page 6 for individual components, including laser components.



A

M3 size resilient: **OL 694 RNX**



B

M3 size resilient: **OL 694 KSNX**



C

M3 size resilient: **OL 694 SNX**



A

M3 size resilient: **OL 694 R**



B

M3 size resilient: **OL 694 KS**



C

M3 size resilient: **OL 694 S**



A

M3 size resilient: **OL 694 RPR**



B

M3 size resilient: **OL 694 KSPR**



C

M3 size resilient: **OL 694 SPR**

CEKA KITS

NEW

*



**



LABORATORY TOOL KIT* **# REVAX LK 3**

Contents: H 1, H 13, H 14, H 30, H 16, H 2, A 1, AF 69, P 7, P 8, 694 C, RE 0031, RE H 1, RE H 10, RE H 13, RE H 14, RE H 16

REDUCED LABORATORY TOOL KIT **# REVAX RK 3**

Contents: H 1, H 13, H 14, H 16, H 2, P 7, P 8, RE H 1, RE H 10, RE H 13, RE H 14, RE H 16, RE H 2, RE H 5, RE P 2/5, RE P 7

DENTIST TOOL KIT **# REVAX OK 2**

Contents: AF 69, A 1, H 30, RE H 5, RE H 79, RE 0031, 694 C, 724 C

SPRING PIN KIT M2 & M3** **# REVAX SP 1**

Contents: RE 0031, RE 0031 L, RE 0031 195, 694 C, 694 C/L, 694 C/195, 694 C/202, 694 C/209, AF 69, A 1, RE H 5, CB1, H 30

SPRING PIN KIT M2 **# REVAX SP 2**

Contents: RE 0031, RE 0031 L, RE 0031 195, A 1, RE H 5, CB1, H 30

SPRING PIN KIT M3 **# REVAX SP 3**

Contents: 694 C, 694 C/L, 694 C/195, 694 C/202, 694 C/209, AF 69, A 1, RE H 5, CB1, H 30

MORE AT: WWW.CEKA-PRECILINE.COM

CEKA INDIVIDUAL COMPONENTS*

① 

M2 size: **RE 0085**
M3 size: **694 AR**
Retention part for bonding with CEKA SITE. Including duplicating dummy.
TITANAX alloy - 2 + 2 pcs

② 

M2 size: **RE 0075**
M3 size: **694 AKS**
Retention part for winged acrylic fixation
TITANAX alloy - 2 pcs

③ 

M2 size: **RE 0095**
M3 size: **694 AKS2**
Retention part for circular acrylic fixation
TITANAX alloy - 2 pcs

④ 

M2 size: **RE 0061**
M3 size: **694 AS**
Retention part for soldering with CEKA SOL
PALLAX alloy - 2 pcs

NEW 

M2 size: **RE 0065**
M3 size: **694 AL**
Retention part for the laser technique
PURE TITANIUM - 2 pcs

⑤ 

M2 size: **RE 0031**
M3 size: **694 C**
Male spring pin. Use CEKA BOND to prevent gradual unthreading.
PALLAX alloy - 2 pcs + 2 space maintainers



M3 size only: **724 C**
Rigid male spring pin. Use CEKA BOND to prevent gradual unthreading.
PALLAX alloy - 2 pcs

⑥ 

M2 size: **RE 0096**
M3 size: **694 B**
Space maintainer of 0.3 mm thickness
STAINLESS STEEL - 10 pcs

⑦ 

M2 size: **RE 0100 TI**
M3 size: **OL 0200 TI**
Female for bonding with CEKA SITE.
Use any dental alloy.
TITANAX alloy - 1 + 3 pcs

⑧ 

M2 size: **RE 0100 NP**
M3 size: **OL 0200 NP**
Female for the cast-on technique with non-precious alloys
NOPRAX alloy - 1 + 3 pcs

⑨ 

M2 size: **RE 0100 IR**
M3 size: **OL 0200 IR**
Female for the cast-on technique with precious alloys
IRAX alloy - 1 + 3 pcs

⑩ 

M2 size: **RA 0061**
M3 size: **691 D**
Base ring for soldering with
CEKA SOL FILIGRAN
PALLAX alloy - 1 pc

⑪ 

M2 size: **RA 0063**
M3 size: **693 D**
Base ring for direct casting with precious alloys
IRAX alloy - 1 pc

NEW 

M2 size: **RA 0064**
M3 size: **693 DNX**
Base ring for direct casting with non-precious alloys
NOPRAX alloy - 1 pc

⑫ 

M2 size: **RA 0091**
M3 size: **691 A**
Female for retention in acrylic overdentures
TITANAX alloy - 1 pc

⑬ 

M2 size: **RA 0055**
M3 size: **691 E**
Large space maintainer for post-copings
TINFOIL - 1 pc

⑭ 

M2 size: **RE 0700 TI**
M3 size: **OL 0800 TI**
Female for bonding with CEKA SITE in post-copings or in bar constructions
TITANAX alloy - 1 + 1 pc

⑮

M3 size only: **OL 694 DNX**
Female for the cast-on technique with non-precious alloys
NOPRAX alloy - 1 pc

⑯ 

M3 size only: **OL 694 D**
Female for the cast-on technique with precious alloys
IRAX alloy - 1 pc

⑰ 

M3 size only: **OL 694 DPR**
Female for the cast-on technique with any dental alloy
PLATIRAX alloy - 1 pc

⑱ 

M3 size only: **334 A**
Retention nut for 334 C
TITANAX alloy - 1 pc

⑲ 

M3 size only: **334 C**
334 male spring pin for the bolt-on technique
PALLAX alloy - 1 pc

⑳ 

M3 size only: **334 F**
Burn-out plastic ring for the bolt-on technique
1 pc

㉑ 

M2 size: **RE 0083 IR**
M3 size: **703 A**
CEKA ID retention part for soldering with CEKA SOL
IRAX alloy - 1 pc

㉒ 

M2 size: **RE 0053 IR**
M3 size: **703 D**
Female for CEKA ID for the cast-on technique with precious alloys
IRAX alloy - 1 pc

* See page 19 for table of alloys and materials.

OVERSIZED MALE SPRING PINS*



M2 size spring pin (PALLAX alloy) with larger head (ø 1.95 mm)
 # RE 0031 195 **2 pcs**



M2 size spring pin (gold-plated PALLAX alloy) extra long (0.3 mm longer)
 # RE 0031 L **2 pcs**



M3 size spring pin (PALLAX alloy) with larger head (ø 1.95 mm, 2.02 mm, 2.09 mm)
 # 694 C/195 **2 pcs**
 # 694 C/202 **2 pcs**
 # 694 C/209 **2 pcs**



M3 size spring pin (PALLAX alloy) extra long (0.3 mm longer)
 # 694 C/L **2 pcs**



Rigid M3 size spring pin (PALLAX alloy) with larger head (ø 1.95 mm, 2.02 mm, 2.09 mm)
 # 724 C/195 **2 pcs**
 # 724 C/202 **2 pcs**
 # 724 C/209 **2 pcs**



M3 size spacer disc (PALLAX alloy), to lengthen the spring pin. Use maximum 5 discs. Use also CEKA BOND.
 # AF 69 0.05 mm thick **10 pcs**

Old spring pins are still available. See separate leaflet no 98.E or www.ceka-preciline.com.

UNIVERSAL ACCESSORIES M2 AND M3 SIZE



M2 & M3 size: **RE H 5**
Laboratory key



M2 & M3 size: **A 1**
Activating tool and screw-driver



M2 & M3 size: **RE H 79**
Instrument for retention reduction



M2 & M3 size: **H 30**
Set of 4 profiles to indicate the diameter of the females

CEKA SITE
CEKA BOND
See page 18.



M2 & M3 size: **P 8**
Paralleling mandrel and sandblasting tool for TITANAX females



M2 & M3 size: **RE P 2/5**
Paralleling mandrel for females

M2 & M3 size: **RE H 20**
Diamond burr for the bonding technique
2 pcs

M2 & M3 size: **RE H 10**
Carbide burr with smooth stop for the metal sleeve

M2 & M3 size: **P 9**
Paralleling mandrel for CEKA ID

ACCESSORIES M2 OR M3 SIZE

M2 size: **RE P 7**
M3 size: **P 7**
Paralleling mandrel for green female keepers

M2 size: **RE P 4**
M3 size: **P 4**
Paralleling mandrel for axial base rings

M2 size: **RE H 2** **2 pcs**
M3 size: **H 2** **2 pcs**
Impression tool for axial attachments, to be used with the (RE) H 14

M2 size: **RE H 14** **2 + 2 pcs**
M3 size: **H 14** **2 + 2 pcs**
Impression tool for females, to be used with the spring pin, the space maintainer and the (RE) H 13.
Combine also with the (RE) H 2.

NEW

M2 size: **RE H 17** **2 pcs**
M3 size: **H 17** **2 pcs**
Protection cap in TITANAX alloy for threading into the base ring

M2 size: **RE H 13** **2 pcs**
M3 size: **H 13** **2 pcs**
Female analogue for the repair and rebasing technique

M2 size: **RE H 1** **2 pcs**
M3 size: **H 1** **2 pcs**
Dummy spring pin with incorporated space maintainer

M2 size: **RE H 16** **2 pcs**
M3 size: **H 16** **2 pcs**
Soldering accessory for extracoronal retention parts

M2 size: **RE H 4** **2 pcs**
M3 size: **H 4** **2 pcs**
Soldering accessory for axial base rings


* See page 19 for table of alloys and materials.

PRECI-VERTIX®

Extracoronal slide attachment with burn-out plastic male.
The elastic females are easily replaceable and available in three friction levels.
PRECI-VERTIX has easy procedures and minimal servicing.
It gives excellent patient comfort and can be used with anterior abutments for optimum esthetics.
To be used with any dental alloy. All tools included. Over 30 years of clinical experience.



PRECI-VERTIX STANDARD
 # 1801 **KIT**
6 males with attached paralleling mandrel
6 yellow females
1 insertion tool

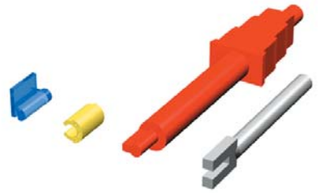

PRECI-VERTIX P
 # 1811 **KIT**
6 angled males for best abutment and tissue contours
6 yellow females
1 paralleling mandrel
1 insertion tool



PRECI-VERTIX P COMBI
 # 1812 **KIT**
2 angled males for best abutment and tissue contours
2 + 2 + 2 females
2 castable housings
1 paralleling mandrel
1 insertion tool

INDIVIDUAL COMPONENTS* & ACCESSORIES



PRECI-VERTIX FEMALE
Polypropylene
Height: 5 mm (may be reduced) - ø 3.30 mm
Available in three friction levels:

normal friction:
 # 1802 yellow **6 pcs**
 # 1802B yellow **50 pcs**
reduced friction:
 # 1805 white **6 pcs**
 # 1805B white **50 pcs**
increased friction:
 # 1806 red **6 pcs**
 # 1806B red **50 pcs**
combination package:
 # 1807 **2 + 2 + 2 pcs**



PRECI-VERTIX STANDARD MALE
Male with extension as paralleling mandrel.
Remove the extension with a hot instrument after positioning. Cast in any hard alloy.
Round the occlusal aspect of the male after processing.
Height: 5.9 mm (may be reduced up to 3 mm)
ø 1.8 mm
 # 1803 **6 pcs**
 # 1803B **50 pcs**



PRECI-VERTIX P MALE
Male with 45° inclination for maximum incorporation possibilities and a periodontal friendly design. Cast in any hard alloy.
Round the occlusal aspect of the male after processing.
Height: 5.9 mm (may be reduced to 3 mm)
ø 1.8 mm
 # 1813 **6 pcs + paralleling mandrel**
 # 1813B **50 pcs**



PRECI-VERTIX PARALLELING MANDREL
Adjustable metal paralleling mandrel for frequent users of PRECI-VERTIX P and AT
Groove width: 1.1 mm
ø paralleling mandrel: 2.35 mm
 # 1845 **1 pc**



PRECI-VERTIX HOUSING
Thin plastic female housing for casting and to incorporate with the soldering, bonding or acrylic fixation technique.
To be adjusted according to the male length.
Process the housing with the white female and replace the latter after finishing with the final female. Sandblast the inside with coarse aluminous oxide.
Height: 7 mm - Length: 12 mm
Extension height: 5 mm - ø 4.1 mm
 # 1808 **6 pcs**
 # 1808B **50 pcs**



PRECI-VERTIX INSERTION TOOL
Double sided insertion tool for the 3 female types
 # 1804 **1 pc**



PRECI-VERTIX MALE ANALOGUE
Brass male analogue, for repairs or for making a new prosthesis
Height: 9 mm - ø 1.8 mm
 # 1819 **2 pcs**



* See page 19 for table of alloys and materials.

PRECI-VERTIX® AT

H 4.4 mm - ø 2.6 mm

Extracoronal friction slide attachment.
The friction of the titanium female can be progressively adjusted with a screw and a plastic friction part.
Good lateral stabilization.
Supplied with a duplicating dummy for bonding with CEKA SITE (female) and a paralleling mandrel.



PRECI-VERTIX AT K
 # 1820 **KIT**
1 burn-out plastic male
1 titanium female with screw and friction part
1 duplicating dummy
1 paralleling mandrel



PRECI-VERTIX AT M
 # 1821 **KIT**
1 cast-on male for precious alloys
1 titanium female with screw and friction part
1 duplicating dummy
1 paralleling mandrel

INDIVIDUAL COMPONENTS* & ACCESSORIES



PRECI-VERTIX AT DUPLICATING DUMMY
Plastic duplicating dummy to prepare the refractory model of the cast metal partial frame for bonding with CEKA SITE
Height: 4.5 mm - Length: 6.1 mm
Width: 3.1 mm
 # 1827 **6 pcs**



PRECI-VERTIX PARALLELING MANDREL
Adjustable metal paralleling mandrel for frequent users of PRECI-VERTIX AT and P
Groove width: 1.1 mm
ø paralleling mandrel: 2.35 mm
 # 1845 **1 pc**




PRECI-VERTIX AT FEMALE
Female in TITANAX alloy with screw and friction part.
May be reduced by 1 mm at tissue side.
Height: 4.4 mm - Length: 6 mm
Width: 2.6 mm
 # 1826 **1 pc**



PRECI-VERTIX AT SCREW-DRIVER
Sterilizable screw-driver for the activating screw
Blade: 0.3 mm - Total length: 30 mm
Handle length: 17.5 mm - ø 1.1 mm
 # 1829 **1 pc**



PRECI-VERTIX AT ACTIVATING SCREW
Screw in TITANAX alloy to activate the friction part
Height: 3 mm - Slot: 0.3 mm - M 1.4
 # 1828 **1 pc**

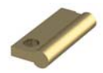

PRECI-VERTIX AT MALE ANALOGUE
Brass male analogue, for repairs or when making a new prosthesis
Height: 9 mm - ø 1.8 mm
 # 1819 **2 pcs**



PRECI-VERTIX AT FRICTION PART
Removable plastic friction part.
Yellow for adjustment of normal friction.
Red for adjustment of increased friction.
Acetal homopolymer - Height: 3.75 mm
 # 1822 yellow **6 pcs**
 # 1822R red **6 pcs**



PRECI-VERTIX AT K MALE
Burn-out plastic male.
May be reduced by 1 mm at tissue side.
Cast in any hard alloy.
Height: 4 mm - ø 1.78 mm
 # 1824 **6 pcs**



PRECI-VERTIX AT M MALE
Male in the high-fusing IRAX alloy for direct casting with precious alloys.
May be reduced by 1 mm at tissue side.
Height: 4 mm - ø 1.78 mm
 # 1823 **1 pc**

CEKA SITE
CEKA BOND
See page 18.


* See page 19 for table of alloys and materials.

MORE AT: WWW.CEKA-PRECILINE.COM

PRECI-SAGIX

Mini: H 3,1 mm - ø 3,1 mm
Standaard: H 4,1 mm - ø 4,2 mm

Sagittal ball attachment with segmented female for partial dentures and implant bar applications.
The replaceable plastic female in three retention levels is incorporated directly into the framework.
The male is in plastic and cast in a hard alloy, in metal and cast-on with a non-precious alloy, or threaded into a cast-on base ring.
PRECI-SAGIX is available in two ball sizes: Mini 1.7 mm and Standard 2.2 mm.

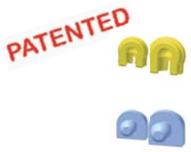





PRECI-SAGIX INTRO KIT P
Introduction kit with **castable plastic males**, the necessary tools and females in three different retention levels for Mini size 1.7 and Standard size 2.2.
Cast the males in any hard dental alloy.

22-piece kit:
2 + 2 castable **plastic** males
2 + 2 yellow females (normal retention)
2 + 2 white females (reduced retention)
2 + 2 red females (increased retention)
2 + 2 duplicating dummies/processing females
1 universal paralleling mandrel
1 universal insertion tool

Intro kit castable male: # 1301

PRECI-SAGIX INTRO KIT M
Introduction kit with **cast-on NOPRAX males**, the necessary tools and females in three different retention levels for Mini size 1.7 and Standard size 2.2.
Cast-on technique for non-precious alloys.

22-piece kit:
2 + 2 cast-on **NOPRAX** males
2 + 2 yellow females (normal retention)
2 + 2 white females (reduced retention)
2 + 2 red females (increased retention)
2 + 2 duplicating dummies/processing females
1 universal paralleling mandrel
1 universal insertion tool

Intro kit cast-on male: # 1303

INDIVIDUAL COMPONENTS* & ACCESSORIES



PRECI-SAGIX FEMALE
Polyacetal

Mini size 1.7: height: 3.1 mm - ø 3.1 mm
yellow - normal retention:
 # 1311 6 pcs # 1311B 50 pcs
white - reduced retention:
 # 1312 6 pcs # 1312B 50 pcs
red - increased retention:
 # 1313 6 pcs # 1313B 50 pcs
combination package:
 # 1314 2 + 2 + 2 pcs

Standard size 2.2: height: 4.1 mm - ø 4.2 mm
yellow - normal retention:
 # 1361 6 pcs # 1361B 50 pcs
white - reduced retention:
 # 1362 6 pcs # 1362B 50 pcs
red - increased retention:
 # 1363 6 pcs # 1363B 50 pcs
combination package:
 # 1364 2 + 2 + 2 pcs

PRECI-SAGIX DUPLICATING DUMMY / PROCESSING FEMALE
Mini size 1.7:
 # 1316 6 pcs # 1316B 50 pcs
Standard size 2.2:
 # 1366 6 pcs # 1366B 50 pcs

PRECI-SAGIX MALE P
Castable spherical plastic male.
Cast in a hard alloy.
Mini size 1.7:
 # 1321 6 pcs # 1321B 50 pcs
Standard size 2.2:
 # 1371 6 pcs # 1371B 50 pcs

PRECI-SAGIX MALE M
NOPRAX male for cast-on technique with non-precious alloys
Mini size 1.7:
 # 1322 2 pcs # 1322B 20 pcs
Standard size 2.2:
 # 1372 2 pcs # 1372B 20 pcs

PRECI-SAGIX MALE V
NOPRAX Standard male ø 2.2 mm with M2 threading and 0.9 mm internal hex.
Use the IMP-XS-042 screw-driver.
 # 1373 1 pc



PRECI-SAGIX BASE RING
For cast-on technique in IRAX (with precious alloys) or in NOPRAX (with non-precious alloys).
M2 threading for the 1373 Standard male
ø 2.2 mm
 # RA 0063 IRAX ALLOY 1 pc
 # RA 0064 NOPRAX ALLOY 1 pc

NEW



PRECI-SAGIX HOUSING
Cr-Co female housing (Co 65 - Cr 28 - Mo 6) for incorporation into acrylic resin, for soldering with CEKA SOL or for the laser technique.
Mini size 1.7:
height: 3.3 mm - width: 3.8 mm - length: 9.1 mm
 # 1325 2 pcs # 1325B 20 pcs
Standard size 2.2:
height: 4.5 mm - width: 4.9 mm - length: 9.5 mm
 # 1375 2 pcs # 1375B 20 pcs



PRECI-SAGIX PARALLELING MANDREL MALE
Universal paralleling mandrel for the PRECI-SAGIX males
 # 1300 P for ø 1.7 and 2.2 1 pc



PRECI-SAGIX SCREW-DRIVER FOR THREADED MALE
Total length: 26 mm - External hex: 0.9 mm
 # IMP-XS-042 1 pc



CEKA BOND
Anaerobic adhesive to prevent gradual unthreading of the male.
 # CB1 5 ml



PRECI-SAGIX INSERTION TOOL
Universal insertion tool for the PRECI-SAGIX females
 # 1323 for ø 1.7 and 2.2 1 pc



PLASTICWAX 8° PROFILE
Height: 4.5 mm - Length: 60 mm
Width top: 1.9 mm / bottom: 2.8 mm
 # IMP-CK-014 4 pcs



PRECI-SAGIX ANALOGUE
Model analogue of male
 # 1331 for ø 1.7 2 pcs
 # 1381 for ø 2.2 2 pcs

* See page 19 for table of alloys and materials.

PRECI-BALL

H 3.55 mm - 4.35 mm
ø 3.5 mm

Resilient adjustable ball attachment. The 2.25 mm ball diameter makes PRECI-BALL compatible with most implant systems.
The male can be soldered, threaded or cast. The female is available in gold alloy or titanium.
See page 12 for individual components.



PRECI-BALL STANDARD OR / OR
 # 1201
Gold female for acrylic fixation and gold male for soldering with CEKA SOL FILIGRAN to the post-coping



PRECI-BALL STANDARD TI / PA
 # 1203
Titanium female for acrylic fixation and palladium male for soldering with CEKA SOL FILIGRAN to the post-coping



PRECI-BALL PC
 # 1205
Titanium female for acrylic fixation and burn-out plastic male for casting



PRECI-BALL CV
 # 1206
Gold female for acrylic fixation and removable titanium male with high-fusing base ring for direct casting with precious alloys

PRECI-CLIX RADICULAR

Chairside ball attachment for acrylic resin overdentures with replaceable plastic female in three retention levels.
The female can be incorporated into acrylic resin or directly into the framework.
The titanium post is available in two versions: with removable male or with the spherical male in one piece.
See page 12 for individual components.



PRECI-CLIX RTI # 1291
Female housing for incorporation into acrylic resin. Cementable titanium post with removable titanium male.



PRECI-CLIX RP # 1292
Replaceable plastic female and duplicating dummy for cast metal partial frames. Cementable titanium post with removable titanium male.

NEW



PRECI-CLIX RCT # 1293
Female housing for incorporation into acrylic resin. Cementable titanium post with titanium male.

NEW



PRECI-CLIX RCP # 1294
Replaceable plastic female and duplicating dummy for cast metal partial frames. Cementable titanium post with titanium male.





NEW



PRECI-CLIX DENTIST KIT
 # 1299
Contains 2 radicular PRECI-CLIX RTI ball attachments, all necessary accessories and spare parts.
Set:

2 x 1291	radicular PRECI-CLIX RTI ball attachments
1 x 1226	set of 3 reamers
1 x 1222	insertion tool
1 x IMP-XS-042	screw-driver
1 x RE H 2	impression tool (2 pcs)
1 x RE H 14	model analogue (2 pcs)
1 x 1234	replacement females (3 x 2 pcs)
1 x CEKA BOND	anaerobic adhesive to prevent gradual unthreading of the male part

PRECI-CLIX RC DENTIST KIT
 # 1297
Contains 2 radicular PRECI-CLIX RCT ball attachments, all necessary accessories and spare parts.
Set:


2 x 1293	radicular PRECI-CLIX RCT ball attachments
1 x 1226	set of 3 reamers
1 x 1222	insertion tool
1 x 1201 D	male analogue
1 x 1234	replacement females (3 x 2 pcs)

* See page 19 for table of alloys and materials.

PRECI-CLIX AXIAL

H 3.4 mm - 4.1 mm
ø 3.55 mm - 4 mm

Axial ball attachment for acrylic resin overdentures with replaceable plastic female in three retention levels.
The female can be incorporated into acrylic resin or directly into the framework.
The male can either be threaded or cast in any hard alloy.

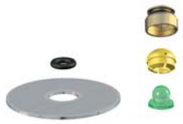

PRECI-CLIX ATI # 1251
Female housing for incorporation into acrylic resin.
Plastic male for incorporation into post-copings.
Cast in any hard alloy.



PRECI-CLIX AP # 1252
Female duplicating dummy for cast metal partial frames.
Plastic male for incorporation into post-copings.
Cast in any hard alloy.



PRECI-CLIX CV # 1261
Female housing for incorporation into acrylic resin.
Removable titanium male with high-fusing base ring for casting with precious alloys.

PRECI-CLIX EXTRACORONAL

H 5.1 mm - 5.3 mm
ø 3.4 mm

Periodontal friendly ball attachment for extracoronal and bar constructions with replaceable plastic female in three retention levels.
The female can be incorporated into acrylic resin or directly into the framework.
The male can be incorporated into the preferred profile by bonding with CEKA SITE, as a threaded component or as a burn-out component.
To be used with any hard alloy. 1 attachment including 4 profiles/package



PRECI-CLIX TI-TI # 1271
Titanium housing for incorporation into acrylic resin.
Titanium male for cementing in the cast male keeper.

NEW



PRECI-CLIX TI-V # 1276
Titanium housing for incorporation into acrylic resin.
Removable titanium male and cementable threading to incorporate.



PRECI-CLIX TI-P # 1281
Titanium housing for incorporation into acrylic resin.
Plastic male to assemble with the male keeper for casting in any hard alloy.



PRECI-CLIX P-TI # 1272
Female duplicating dummy for cast metal partial frames.
Titanium male for cementing in the cast male keeper.

NEW



PRECI-CLIX P-V # 1277
Female duplicating dummy for cast metal partial frames.
Removable titanium male and cementable threading to incorporate.



PRECI-CLIX P-P # 1282
Female duplicating dummy for cast metal partial frames.
Plastic male to assemble with the male keeper for casting in any hard alloy.

PRECI-BALL/CLIX - INDIVIDUAL COMPONENTS* & ACCESSORIES



PRECI-BALL FEMALE
Female in TITANAX or in ORAX alloy for acrylic fixation
Height: 2.8 mm - ø 3 mm

# 1205 A	TITANAX ALLOY	1 pc
# 1205 AB	TITANAX ALLOY	20 pcs
# 1201 A	ORAX ALLOY	1 pc
# 1201 AB	ORAX ALLOY	20 pcs
# 1201 B	processing ring	1 pc



PRECI-BALL/CLIX LARGE SPACE MAINTAINER
TINFOIL - Thickness: 0.4 mm - ø 14 mm
1201 E 1 pc



PRECI-CLIX HOUSING
Female housing in an anodized titanium alloy for fixation in acrylic resin dentures
Height: 2.65 mm - ø 4 mm

# 1235	6 pcs
# 1235B	50 pcs



PRECI-CLIX FEMALE
Polyacetal - Height: 2.65 mm - ø 3.55 mm
Normal retention:

# 1231	yellow	6 pcs
# 1231B	yellow	50 pcs

Reduced retention:

# 1232	white	6 pcs
# 1232B	white	50 pcs

Increased retention:

# 1233	red	6 pcs
# 1233B	red	50 pcs

Combination package:
1234 2 + 2 + 2 pcs

* See page 19 for table of alloys and materials.



PRECI-CLIX DUPLICATING DUMMY
Plastic duplicating dummy for cast metal partial frames. No need for a space maintainer!
Height: 2.8 mm - ø 3.55 mm
1236 6 pcs



PRECI-CLIX SPACE MAINTAINER
Elastic space maintainer, to be placed between male and female during processing
Rubber - Height: 0.7 mm - ø 1.8 mm
1251 B 6 pcs

NEW



PRECI-CLIX THREADS TI
M2 threading in TITANAX alloy for cementing with CEKA SITE.
Combine with the 1206 C PRECI-CLIX male V.
Height: 1.9 mm - ø 2.3 mm
RE 4600 TI 1 pc



PRECI-CLIX MALE KEEPER
4 burn-out plastic male keepers for extracoronal and bar constructions
Adaptable in form and height
Internal ø 2.35 mm - External ø 3.4 mm
RE 04 TI 1 x 0°, 1 x 30°, 1 x 45°, 1 x 60°

NEW



PRECI-CLIX POST
Cementable post in TITANAX alloy with threaded base for the 1206 C PRECI-CLIX male V or with attached ball
Post length: 7 mm (may be reduced) - ø 1.35 mm
ø base: 3.4 mm - Threads: M2 - ø ball: 2.25 mm
1291 C with threads 1 pc
1293 C with ball 1 pc



PRECI-CLIX-R BURR SET
Set of 3 calibrated burrs for right angle.
Used for root preparation prior to cementation of the PRECI-CLIX post.

# 1226	Set: 1 predrilling burr, 1 cavity burr, 1 precision reamer
# 1227	1 predrilling burr
# 1228	1 cavity burr
# 1229	1 precision reamer



PRECI-BALL MALE
Spherical male in ORAX alloy (yellow) or PALLAX alloy (white) for soldering with CEKA SOL FILIGRAN
Height: 2.43 mm – ø 2.25 mm – ø base: 3.2 mm
1201 C ORAX ALLOY 1 pc
1203 C PALLAX ALLOY 1 pc



PRECI-BALL/CLIX MALE PA
Spherical burn-out plastic male for incorporation into post-copings and for casting
Height: 2.4 mm - ø 2.25 mm - ø base: 3.2 mm
1205 C 1 pc
1205 CB 50 pcs

NEW



PRECI-BALL/CLIX MALE V AND V LONG
Spherical, removable male in TITANAX alloy for incorporation into an M2 base ring on post-copings
Height male V: 3.05 mm, male V long: 3.35 mm
ø 2.25 mm - Threads: M2 - Internal hex: 0.9 mm
1206 C 1 pc
1206 CL 0.3 mm extra long 1 pc



PRECI-BALL/CLIX BASE RING
M2 base ring in the high-fusing IRAX or NOPRAX alloy to incorporate the 1206 C PRECI-CLIX male V. For direct casting or for soldering.
Height: 1.3 mm - ø 3.4 mm - Threads: M2
RA 0063 IRAX ALLOY 1 pc
RA 0064 NOPRAX ALLOY 1 pc



PRECI-CLIX MALE PE
Spherical plastic male for incorporation into the male keeper and for casting
Height: 3.95 mm - ø 2.25 mm
1281 C 1 pc



PRECI-CLIX MALE TI
Spherical titanium male for cementation into a cast male keeper. Use CEKA SITE.
Height: 3.95 mm - ø 2.25 mm
1271 C 1 pc



IMPRESSION TOOL M2
Impression and cementation tool for the post/base ring. Total height: 7 mm - Threads: M2
RE H 2 2 pcs



MODEL ANALOGUE M2
Model analogue for the post/base ring.
To be used with the RE H 2 and the male.
Total height: 6 mm - Threads: M2
RE H 14 2 pcs



PRECI-BALL/CLIX CUP BURR
Burr to finish the 1205 C PRECI-BALL/CLIX male PA and the 1281 C PRECI-CLIX male PE
ø 2.35 mm
1230 2 pcs



PRECI-BALL/CLIX PARALLELING MANDREL MALE
Adjustable paralleling mandrel for males
ø 2.35 mm
1201 P 1 pc



PRECI-BALL/CLIX PARALLELING MANDREL FEMALE
Adjustable paralleling mandrel for females.
To be used with diverging males.
ø 2.35 mm
1211 P 1 pc



PRECI-BALL/CLIX PARALLELING MANDREL BASE RING
Paralleling mandrel for M2 base rings
ø 2.35 mm - Threads: M2
RE P 4 1 pc



PRECI-CLIX PARALLELING MANDREL MALE KEEPER
Paralleling mandrel for the RE 04 TI male keeper
ø 2.35 mm
RE P 7 1 pc



PRECI-BALL/CLIX SCREW-DRIVER FOR THREADED MALES
External hex: 0.9 mm - Length: 26 mm
IMP-XS-042 1 pc



PRECI-BALL/CLIX MALE ANALOGUE
Male analogue for the working model, including space maintainer
1201 D 2 pcs + 2 space maintainers



PRECI-BALL ACTIVATING TOOL
Activating tool to tighten the PRECI-BALL segments
1220 1 pc



PRECI-BALL DEACTIVATING TOOL
Deactivating tool to spread the PRECI-BALL segments
1221 1 pc



PRECI-CLIX INSERTION TOOL
Insertion tool for the plastic females
Total height: 67.5 mm
1222 1 pc

NEW




PRECI-BALL/CLIX PROTECTION CAP
Protection cap in TITANAX alloy for threading into the base ring. ø 3 mm - Threads: M2
RE H 17 2 pcs

CEKA SITE
CEKA BOND
See page 18.


* See page 19 for table of alloys and materials.

PRECI-HORIX®


Burn-out plastic bar with elastic plastic riders.
The bar attachment allows a rotational and a slight vertical movement by compression of the plastic riders.
The riders are easily replaceable and available in three retention levels.

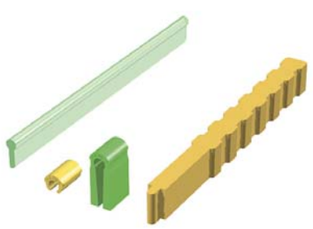

PRECI-HORIX STANDARD
 # 1701 **KIT**
2 plastic bar patterns
6 yellow riders
6 processing riders
1 insertion tool

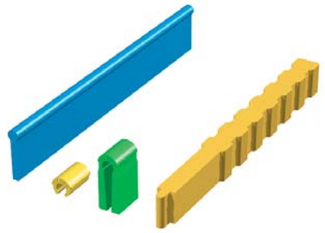

PRECI-HORIX PLASTICWAX
 # 1731 **KIT**
2 PLASTICWAX bar patterns
6 yellow riders
6 processing riders
1 insertion tool

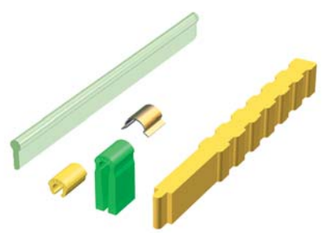

PRECI-HORIX COMBI
 # 1741 **KIT**
1 plastic bar pattern
2 yellow riders
2 processing riders
2 stainless steel housings
1 insertion tool

INDIVIDUAL COMPONENTS* & ACCESSORIES



PRECI-HORIX METAL HOUSING
Ultrathin stainless steel housing with titanium/nitrate layer (gold-coloured) for a perfect fit of the riders in the prosthesis. The housing is polymerized into acrylic resin (with the processing rider).
Do not solder!
Length: 5 mm - Width: 5 mm
Thickness: 0.2 mm
 # 1706 **6 pcs**
 # 1706B **50 pcs**



PRECI-HORIX RIDER
Polypropylene
Height: 2.87 mm - Length: 5 mm
ø 3.3 mm
Available in three retention levels:

normal retention:
 # 1802 yellow **6 pcs**
 # 1802B yellow **50 pcs**
reduced retention:
 # 1805 white **6 pcs**
 # 1805B white **50 pcs**
increased retention:
 # 1806 red **6 pcs**
 # 1806B red **50 pcs**
combination package:
 # 1807 **2 + 2 + 2 pcs**



PRECI-HORIX PROCESSING RIDER
Adapt the plastic processing rider to the bar prior to polymerization.
Length: 5 mm (may be reduced to fit the bar height)
 # 1703 **6 pcs**
 # 1703B **50 pcs**

REDUCED HEIGHT



PRECI-HORIX PLASTIC PATTERN
Burn-out plastic bar pattern which may be reduced to facilitate incorporation into close bite situations. Cast in any hard alloy.
Bar height: 4 mm - Length: 50 mm
ø 1.8 mm
 # 1704 **6 pcs**
 # 1704B **50 pcs**



PRECI-HORIX PLASTICWAX PATTERN
Bar pattern in PLASTICWAX which may be reduced to facilitate incorporation into close bite situations.
Ideal for use with fast set investments. Cast in any hard alloy.
Bar height: 8.25 mm (may be reduced to 3 mm) - Length: 50 mm - ø 1.8 mm
 # 1734 **4 pcs**
 # 1734B **50 pcs**



PRECI-HORIX PATTERN ANALOGUE
Plastic pattern analogue for the working model, to be used in case of repair and relining, or when making a new prosthesis.
Height: 8.25 mm - Length: 50 mm
ø 1.8 mm
 # 1707 **2 pcs**



PRECI-HORIX INSERTION TOOL
Insertion tool for the 3 rider types
 # 1705 **2 pcs**



PRECI-HORIX PARALLELING MANDREL
Paralleling mandrel for bars. Ideal for paralleling multiple pattern sections.
ø 2.35 mm
 # 1708 **1 pc**

* See page 19 for table of alloys and materials.

PRECI-BAR

Bar attachment in two sizes: MINI and STANDARD.
The bar is either rigid (U-shaped) or resilient (oval). A choice of three alloys for the rider.
The bars in gold alloy/palladium alloy can be soldered with CEKA SOL FILIGRAN (yellow) or CEKA SOL W (white), or cast (burn-out plastic).
The metal bars are available in two lengths: 50 mm and 200 mm.

MINI PROFILES

Height: 2.3 mm
Width: 1.6 mm



ORAX* RESILIENT OVAL MINI BAR
 # 1103/B/OR **50 mm**
 # 1103/BL/OR **200 mm**

NEW

PALLAX* RESILIENT OVAL MINI BAR
 # 1103/B/PA **50 mm**
 # 1103/BL/PA **200 mm**

PLASTIC RESILIENT OVAL MINI BAR
 # 1012/B **4 x 50 mm**

MINI RIDERS

Height: 2.8 mm
Width: 4 mm



ORAX* MINI RIDER
 # 1102/H/MR/OR **50 mm**

NEW

PALLAX* MINI RIDER
 # 1102/H/MR/PA **50 mm**

STAINLESS STEEL* MINI RIDER
 # 1102/H/MR/IN **50 mm**

STANDARD PROFILES

Height: 3 mm
Width: 2.2 mm



ORAX* RESILIENT OVAL STANDARD BAR
 # 1101/B/OR **50 mm**
 # 1101/BL/OR **200 mm**

NEW

PALLAX* RESILIENT OVAL STANDARD BAR
 # 1101/B/PA **50 mm**
 # 1101/BL/PA **200 mm**

PLASTIC RESILIENT OVAL STANDARD BAR
 # 1013/B **4 x 50 mm**



ORAX* RIGID U-SHAPED STANDARD BAR
 # 1100/B/OR **50 mm**
 # 1100/BL/OR **200 mm**

NEW

PALLAX* RIGID U-SHAPED STANDARD BAR
 # 1100/B/PA **50 mm**
 # 1100/BL/PA **200 mm**

PLASTIC RIGID U-SHAPED STANDARD BAR
 # 1011/B **4 x 50 mm**

STANDARD RIDERS

Height: 3.5 mm
Width: 5 mm



ORAX* STANDARD RIDER
 # 1100/H/MR/OR **50 mm**

NEW

PALLAX* STANDARD RIDER
 # 1100/H/MR/PA **50 mm**

STAINLESS STEEL* STANDARD RIDER
 # 1100/H/MR/IN **50 mm**

*TABLE OF ALLOYS

ORAX - Yellow - Au 70 - Pt 5 - Ag 12 - Cu 13
Heat treatment: 60 min at 400 °C (752 °F)
Melting range: 910-955 °C (1670-1823 °F)
PALLAX - White - Au 2 - Ag 37 - Pt 9.5 - Pd 37 - Cu 12.5 - Co 2
Melting range: 1055-1130 °C (1931-2066 °F)
STAINLESS STEEL - White - Fe 73 - Cr 18 - Ni 8 - Mo, Cu, C, Mn and Si 1
Do not heat!



PARALLELING MANDRELS
 # 1100 P
Paralleling mandrel for STANDARD bars
 # 1102 P
Paralleling mandrel for MINI bars

PRECI-CLIP

Round bar attachment with adjustable rider.
The rider is available with occlusal or lateral retention for acrylic fixation.
The riders are available in gold alloy or stainless steel.
The profile can be soldered with CEKA SOL FILIGRAN, laser-welded or cast in any hard alloy.



PRECI-CLIP RIDER LA
Adjustable rider with lateral retention for acrylic resin for minimum incorporation possibilities.
Do not solder the rider! Height: 3.65 mm - Length: 3.5 mm - Space maintainer: 0.7 mm
Available in:

# 1105	**ORAX* ALLOY**	**2 pcs + 2 space maintainers**
# 1105AB	**ORAX* ALLOY**	**20 pcs + 20 space maintainers**
# 1106	**STAINLESS STEEL***	**2 pcs + 2 space maintainers**
# 1106AB	**STAINLESS STEEL***	**20 pcs + 20 space maintainers**

NEW



PRECI-CLIP RIDER OC
Adjustable rider with occlusal retention for acrylic resin for minimum incorporation possibilities.
Do not solder the rider! Height: 3.65 mm - Length: 3.5 mm - Space maintainer: 0.7 mm
Available in:

# 1125	**ORAX* ALLOY**	**2 pcs + 2 space maintainers**
# 1125AB	**ORAX* ALLOY**	**20 pcs + 20 space maintainers**
# 1126	**STAINLESS STEEL***	**2 pcs + 2 space maintainers**
# 1126AB	**STAINLESS STEEL***	**20 pcs + 20 space maintainers**



PRECI-CLIP BAR PROFILE
Round bar profile in ORAX* alloy for soldering with CEKA SOL FILIGRAN or in pure titanium for the laser technique
Length: 50 mm - ø 1.8 mm

# 1105/B	**ORAX* ALLOY**	**1 pc**
# 1109/B	**PURE TITANIUM***	**1 pc**



PRECI-CLIP BAR PATTERN
Round plastic bar pattern for casting in any hard alloy between natural abutments and between implants
Length: 50 mm - ø 1.8 mm

# 1014/B	**4 pcs**

PRECI-CLIP PATTERN ANALOGUE
Pattern analogue in hard isolating plastic for the working model, to be used in case of repair and relining, or when making a new prosthesis.
Height: 8.25 mm - Length: 50 mm - ø 1.8 mm

# 1707	**2 pcs**

PRECI-TUBIX / MORTIX

Intracoronal, plastic slide attachment for sectioned bridges and diverging abutments



PRECI-MORTIX N
Conical
Height: 7 mm (may be reduced)
Total width: 2.85 mm
Cone: 2°23'
 # 10002 **6 pcs**



PRECI-TUBIX N
Cylindrical, small
Ideal for use in the anterior region.
Height: 7 mm (may be reduced)
Total width: 2 mm
ø male: 1.05 mm
 # 10001 **6 pcs**



PRECI-TUBIX M
Cylindrical, large, **reversible**
Height: 6 mm (may be reduced)
Total width: 3 mm
ø male: 1.45 mm
 # 10004 **6 pcs**

PRECI-TUBIX MORTIX COMBI
Combination package

# 10005	**2 x 10001, 2 x 10002, 2 x 10004**

* See page 19 for table of alloys and materials.

PRECI-POST M

Rational system with a conicity of 2°30' over the entire useful length.
Form and dimensions according to Dr Mooser. For post-copings and temporary post cores with one or more roots.
Yellow and blue colour codification.



PRECI-POST M SET
 # 2001
1 yellow reamer
1 blue reamer
25 yellow opaque impression posts
25 blue opaque impression posts
25 yellow transparent laboratory posts
25 blue transparent laboratory posts



PRECI-POST M LABORATORY POST
Transparent plastic laboratory posts, to be used also for direct method in the mouth

# 2003GL	**50 pcs**	yellow	tip ø 0.80 mm, cone length 10.70 mm, total length 18.60 mm
# 2003BL	**50 pcs**	blue	tip ø 1.00 mm, cone length 12.30 mm, total length 22.90 mm
# 2003L	**2 x 50 pcs**	combination package yellow & blue	



PRECI-POST M IMPRESSION POST
Opaque plastic impression posts (with retention), to be used with any impression technique (including silverplating and epoxy)

# 2002G	**50 pcs**	yellow	tip ø 0.80 mm, length 14.10 mm
# 2002B	**50 pcs**	blue	tip ø 1.00 mm, length 16.70 mm
# 2002	**2 x 50 pcs**	combination package yellow & blue	



PRECI-POST M AUXILIARY POST
Stainless steel posts (with flat head) for temporary restorations

# 2500	**20 pcs**	yellow	tip ø 0.85 mm, cone length 7.40 mm, total length 11 mm
# 2501	**20 pcs**		tip ø 0.85 mm, cone length 8.50 mm, total length 13 mm
# 2501B	**20 pcs**		tip ø 0.85 mm, cone length 10.70 mm, total length 15 mm
# 2502	**20 pcs**	blue	tip ø 1.05 mm, cone length 10.70 mm, total length 16 mm
# 2502B	**20 pcs**		tip ø 1.05 mm, cone length 11.80 mm, total length 17 mm
# 2503	**20 pcs**		tip ø 1.05 mm, cone length 12.80 mm, total length 19 mm



PRECI-POST M REAMER
Reamer according to the Mooser system

# 2004G	**1 pc**	yellow	tip ø 0.80 mm
# 2004B	**1 pc**	blue	tip ø 1.00 mm

PRECI-POST CC

Cylindro-conical system: 5°30' at the base and then cylindrical over the entire useful length for optimum retention.
For post-copings and temporary post cores with one or more roots.
Yellow, red, blue, and green colour codification.



PRECI-POST CC SET

# 2050	50 yellow posts, 50 red posts, 50 blue posts, 50 green posts
# 2060	20 yellow posts, 15 red posts, 10 blue posts, 5 green posts



PRECI-POST CC PLASTIC POST
Plastic posts for the laboratory

# 2051A	**50 pcs**	yellow	ø 1.40 mm, tip ø 0.65 mm, cone L 4.60 mm, total L 16.60 mm
# 2052B	**50 pcs**	red	ø 1.60 mm, tip ø 0.70 mm, cone L 5.30 mm, total L 18.60 mm
# 2053C	**50 pcs**	blue	ø 1.80 mm, tip ø 0.75 mm, cone L 6.30 mm, total L 20.00 mm
# 2054D	**50 pcs**	green	ø 2.00 mm, tip ø 0.80 mm, cone L 7.05 mm, total L 22.40 mm



PRECI-POST CC AUXILIARY POST
Stainless steel posts for temporary restorations

# 2551A	**20 pcs**	ø 1.40 mm, tip ø 0.70 mm, cone length 4.40 mm, total length 11 mm
# 2552B	**20 pcs**	ø 1.60 mm, tip ø 0.75 mm, cone length 5.30 mm, total length 13 mm
# 2553C	**20 pcs**	ø 1.80 mm, tip ø 0.80 mm, cone length 6.20 mm, total length 15 mm
# 2554D	**20 pcs**	ø 2.00 mm, tip ø 0.80 mm, cone length 7.10 mm, total length 17 mm



PRECI-POST CC REAMER
Cylindro-conical reamer

# 2951A	**1 pc**	yellow	tip ø 0.65 mm
# 2952B	**1 pc**	red	tip ø 0.70 mm
# 2953C	**1 pc**	blue	tip ø 0.75 mm
# 2954D	**1 pc**	green	tip ø 0.80 mm

Ceka Site



Ceka Site is an anaerobic bonding composite.

Mix for 30 seconds (1:1) and apply to sandblasted surfaces within 1 minute.
Hold for 10 minutes.
Excess material can be easily removed.
Store refrigerated but use at room temperature.

 # Ceka Site **2 x 2 g**



Ceka Bond

NEW




Ceka Bond is an anaerobic adhesive to prevent gradual unthreading of prosthetic components.

The components can still be unthreaded by using the proper instruments.

Setting time in the mouth: 3 hours

 # CB1 **5 ml**

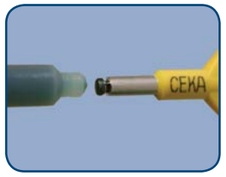

3C-Bond



Bonding agent for porcelain firing to titanium-free, non-precious alloys (Cr-Co-Mo and Ni-Cr).

Avoids greying/greening of porcelain.
Apply the material as a wash opaque.

 # 140-3CB-20 **20 g**



Perma-Ret



Mechanical retention for acrylic and composite teeth.
Ideal for single tooth repairs or in deep bites.

# 1601	**Set:**	**1 wire holder, 3 twist drills** **100 cm of threaded silver wire**
# 1602		**13 x 80 mm of threaded silver wire**
# 1604		**1 wire holder**
# 1606		**3 twist drills**




Plasticwax



Wax with plastic filler.
Ideal for build-up and milling technique.
Available as build-up wax or in 8° profiles (height: 4.5 mm - length: 60 mm - width top: 1.9 mm / bottom: 2.8 mm).

# IMP-CK-045	**box**	**45 g**
# IMP-CK-014	**8° profiles**	**4 pcs**

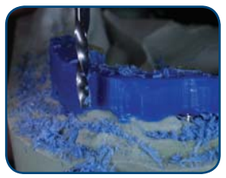

Expando



Expando is a retraction device to protect the lips, cheek, and tongue during dental treatment procedures. Expando helps to relax the patient and also enables a better positioning of the suction device. After placing the expando, the dentist has better access to the treatment area and can keep both hands free to operate. Restorative work can be done in the molar area without chairside assistance. Can be sterilized.

Available in 3 sizes:

# 7301	**small**	**1 pc**
# 7302	**medium**	**1 pc**
# 7303	**large**	**1 pc**



Preci-Profile



Laser welding titanium profiles for wraparound implant bridges
To be used with the CS31 titanium laser welding alloy.
Fast, versatile and economical. Pure titanium grade 2
Height: 2.9 mm - Width: 2 mm - Length: 50 mm - Vickers hardness: 210 ± 20

Straight titanium profile for laser welding
 # IMP-CK-015 **1 pc**
 # IMP-CK-025 **25 pcs**

Curved titanium profile for laser welding
 # IMP-CK-016 **1 pc**
 # IMP-CK-026 **25 pcs**
Curve radius: 20 mm - Heat-treated



Ceka Sol

Ceka Sol is a complete range of low-fusing tube solders for multiple applications.
The incorporated flux provides for easy soldering of dissimilar metals.
Also available as addition solder for laser welding.



Low-fusing tube solder with flux incorporated

# CS00	**Ceka Sol**	**yellow**	**780-820 °C (1436-1508 °F)**	
	ø 1.15 mm 2 x 25 cm		Au 80 - Cu 13 - In 5 - Zn 2	
# CS00 F	**Ceka Sol Filigran**	**yellow**	**780-820 °C (1436-1508 °F)**	
	ø 0.9 mm 2 x 25 cm		Au 80 - Cu 13 - In 5 - Zn 2 **(contains less flux)**	
# CS22	**Ceka Sol W**	**white**	**865-930 °C (1589-1706 °F)**	
	ø 1.15 mm 2 x 25 cm		Au 76.5 - Cu 0.4 - Zn 6.5 - Ni 16.6	
# CS21	**Ceka Sol L**	**yellow**	**615-630 °C (1139-1166 °F)**	
	ø 1.15 mm 2 x 25 cm		Au 42 - Ag 38.5 - Zn 11 - In 8.5	
# CS20	**Ceka Sol G**	**yellow**	**755-820 °C (1391-1508 °F)**	
	ø 1.15 mm 2 x 25 cm		Au 68 - Ag 13.3 - Cu 12.8 - Sn 1 - Zn 4.9	
# CS23	**Ceka Sol CR**	**white**	**1000-1030 °C (1832-1886 °F)**	
	ø 1.15 mm 4 x 25 cm		Au 10 - Ag 37 - Pd 37 - Cu 12 - Zn 4	

PRECIOUS	NON-PRECIOUS	CR-CO	STAINLESS STEEL
X	X	X	X
X	X	X	X
X	X	X	X
X		X	X
X	X		
		X	

Addition alloy for laser welding NEW

# CS31	**Ceka Sol T1**	**white**	**pure titanium**
	ø 0.4 mm 8 x 25 cm		
# CS32	**Ceka Sol OR**	**yellow**	**910-995 °C (1670-1823 °F)**
	ø 0.4 mm 2 x 25 cm		Au 70 - Pt 5 - Ag 12 - Cu 13
# CS33	**Ceka Sol PA**	**white**	**1055-1130 °C (1931-2066 °F)**
	ø 0.4 mm 2 x 25 cm		Au 2 - Pt 9.5 - Ag 37 - Pd 37 - Cu 12.5 - Co 2

Table of Alloys and Materials

Pallax 1055-1130 °C 1931-2066 °F	**For soldering to precious and non-precious alloys** White - Au 2 - Ag 37 - Pt 9.5 - Pd 37 - Cu 12.5 - Co 2
Irax 1400-1460 °C 2552-2660 °F	**For the cast-on technique with precious alloys only** White - Au 60 - Pt 24 - Pd 15 - Ir 1
Noprax 1355-1450 °C 2471-2642 °F	**For the cast-on technique with non-precious alloys only** White - Ni 72 - Cr 17 - Fe 10 - Mn 1
Titanax Do not heat!	**Only for bonding and acrylic fixation** White - Ti 90 - Al 6 - V 4
Orax 910-955 °C 1670-1823 °F	**Heat treatment: 60 min at 400 °C (752 °F)** Yellow - Au 70 - Pt 5 - Ag 12 - Cu 13

Platirax 1790-1820 °C 3254-3308 °F	**For the cast-on technique with all dental alloys** White - Pt 85 - Ir 15
Stainless steel Do not heat!	**Only for bonding and acrylic fixation** White - Fe 66 - Cr 17 - Ni 11 - Mo 2 - other 4
Polypropylene	**Females of Preci-Horix/Vertix** **# 1802, 1805, 1806**
Polyacetal	**Females of Preci-Clix/Sagix** **# 1231, 1232, 1233, 1311, 1312, 1313, 1361, 1362, 1363**
Thermoplastic Elastomer	**Processing ring of Preci-Ball** **# 1201 B**

More at: www.ceka-preciline.com